


EXACTIVE VITAL

BLOOD GLUCOSE METER



QUICK GUIDE

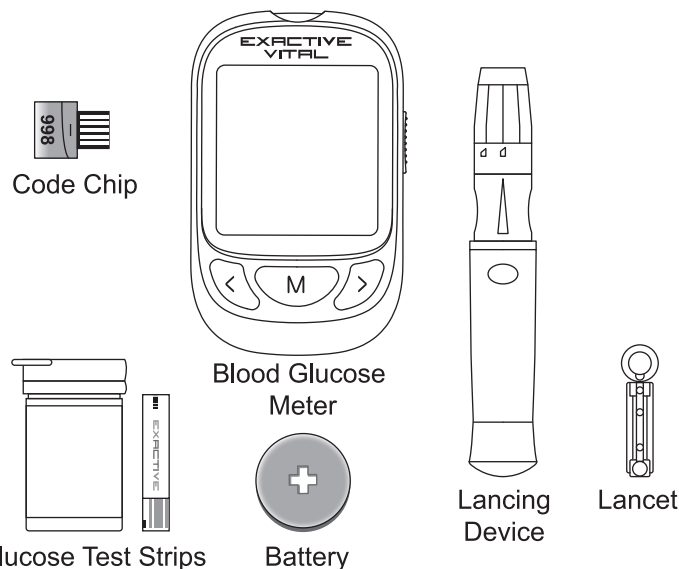

Note: Depending on which Exactive product you purchased, some of the components may need to be purchased separately. Please check the list of components printed on the box for details on which components are included with your purchase.



MicroTech Medical (Hangzhou) Co., Ltd.


Address: No.108 Liuze St., Cangqian, Yuhang District,
Hangzhou, 311121 Zhejiang, P.R.China
Website: <http://www.microtechmd.com>




Lotus NL B.V.
Koningin Julianaplein 10, 1e Verd,
2595AA, The Hague, Netherlands.

1004-FMTL-0201701
Revision 7.0


1 ▶ Insert the battery.

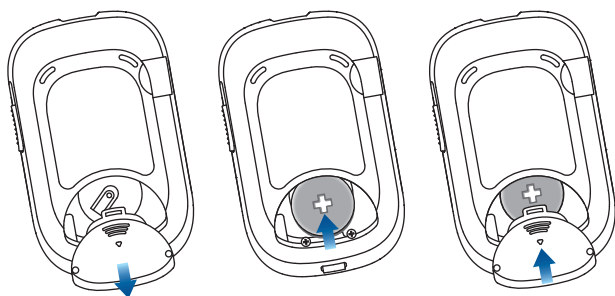


2 ▶ Insert the code chip and compare the number to the number printed on the package, confirm the two are the same.

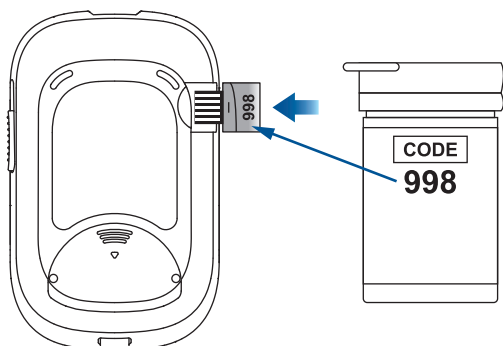

3 ▶ Wash and dry your hands.

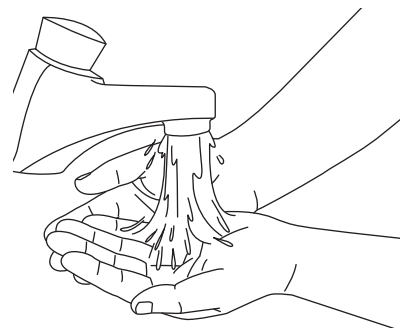


4 ▶ Remove the lancing device cap.

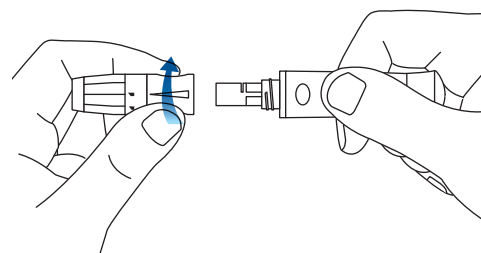


Continued...


- 5** ▶ Insert a lancet into the lancing device.

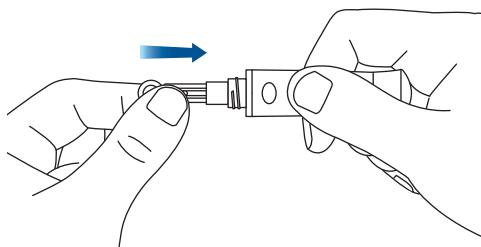


- 6** ▶ Twist off the safety cap and put aside. Do not discard, it will be used later.

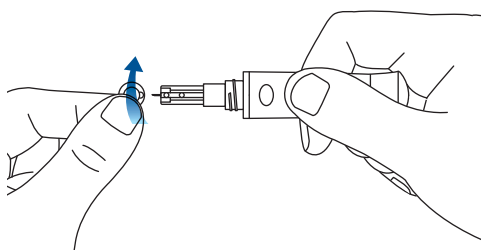


- 7** ▶ Replace the lancing device cap.

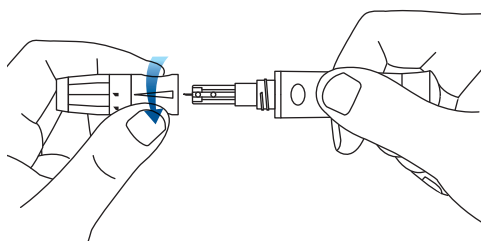


- 8** ▶ Adjust the puncture depth.

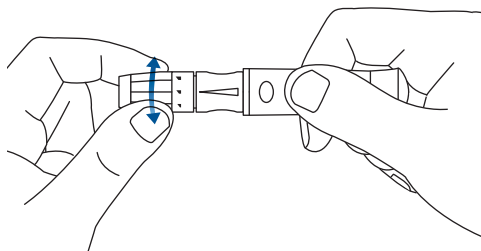


- 9** ▶ Pull back the cocking barrel until you hear a click.

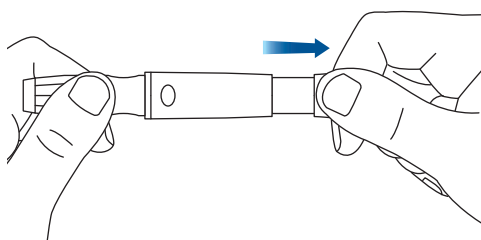


- 10** ▶ Insert a test strip and confirm the code chip number.

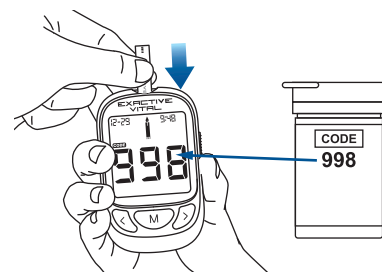


- 11** ▶ Holding the lancing device against the side of the finger, press the release button.

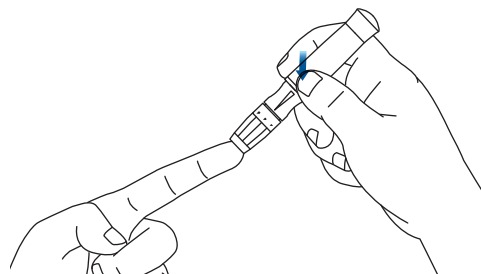


- 12** ▶ Massage the finger to increase the sample size. Wipe away the first drop, and use the second drop to test.

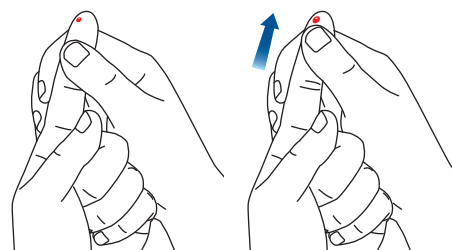

- 13** ▶ Apply blood to the sample area of the test strip. When blood fills the test window, the meter will count down.

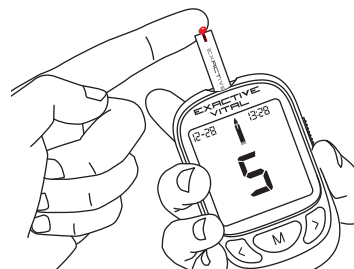

- 14** ▶ After 5 seconds, the result will display.

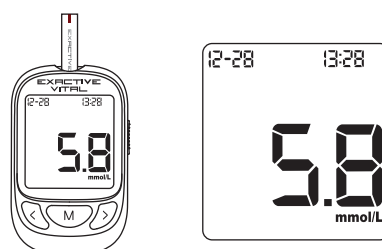

After the test, eject the strip using the test strip ejector, replace the lancet safety cap, and discard the lancet.